User Manual

Model: DM12


Contents

10.1-inch Digital Microscope (EN).....	01
10,1-Zoll-Digitalmikroskop (DE).....	07
Microscope numérique de 10,1 pouces (FR).....	13
Microscopio digitale da 10,1 pollici (IT).....	19
Microscopio digital de 10,1 pulgadas (ES).....	25
10.1インチデジタルマイクロスコープ (JA).....	31
10,1-tommers digital mikroskop (DA).....	37
10,1-tums digital mikroskop (SE).....	43
Digitální mikroskop 10,1 palců (CZ).....	49

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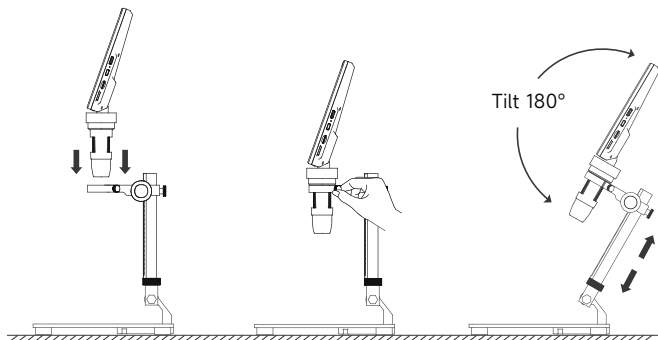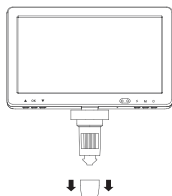
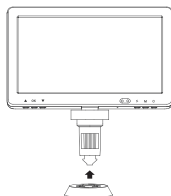
Step2:

Replace the light shield


Step3: When installing, rotate left and right so that the buckle is completely stuck.



Product Description

Parts and functions

Notice: The picture is for reference only, please refer to the actual object

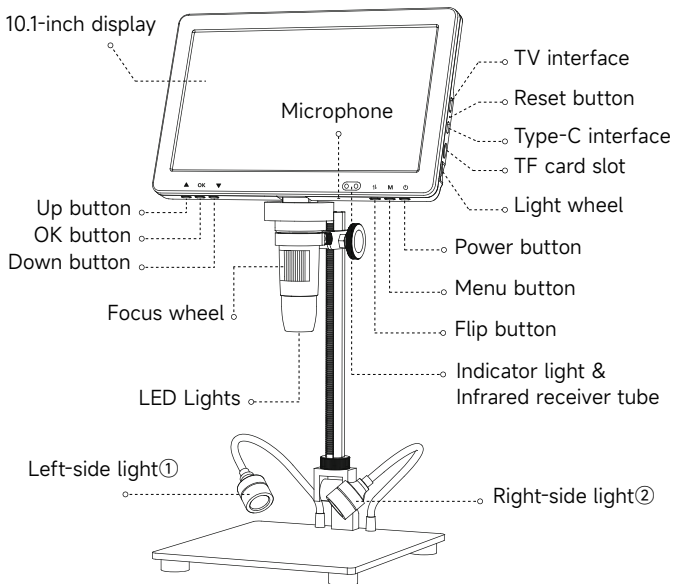Turn on auxiliary lights

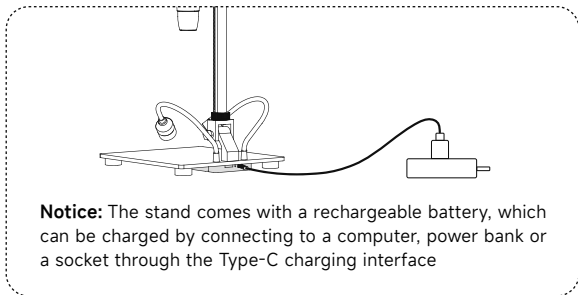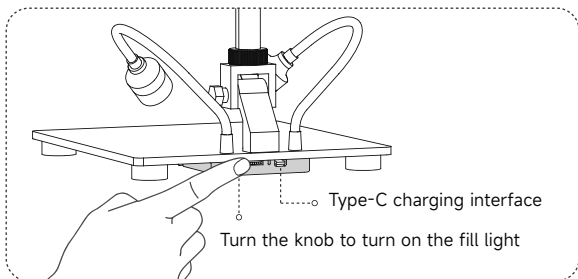Work with PC


Software download

Note: The following 3 methods can be downloaded, please choose the appropriate version according to your computer system.

- ① Please download the applicable software from the following website: **<http://soft.hvscam.com>**
(Suitable for: Windows 7, Windows 8.1, Windows 10, Windows 11 and Mac)
- ② For Windows, search for “HiView” in “Microsoft Store”(“”)
(Suitable for: Windows 8, Windows 10, Windows 11 and above), and you can also be viewed in Camera(“”).
- ③ For Mac OS x 12 or above, search for “HiView” in “App Store”, and you can also be viewed in Photo Booth(“”)

Operating steps

Step 1. Connect the microscope to the computer via Type-C cable.

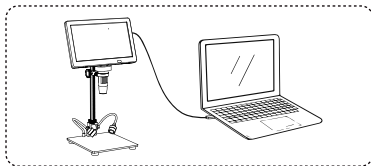
Step 2. Long press the microscope power button to turn on, and the computer side of the successful connection screen pop-up window: Mass Storage/PC Camera/Charge Mode three modes, please select the "PC Camera" mode.

Step 3. Run Hiview, you can watch the image, at this time the photo/video files will be automatically saved in the computer disk. (Note: When connecting to a computer to watch images, please make sure you are using a cable with data transfer function)

Function Introduction of USB Mode

When the device is connected to the computer through Type-C cable, the device can be used in Mass Storage / PC Camera / Charge mode.

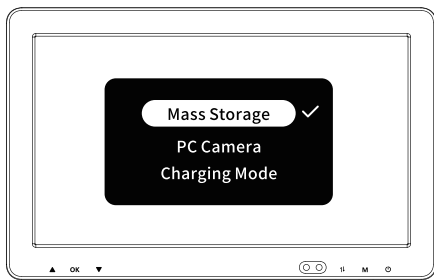①**Mass Storage:** Enter the microscope memory from the computer side, you can view and process the current picture and video files in the TF card, at this time the microscope only Power button and Light wheel available.


②**PC Camera:** After selecting this mode, you can get a larger magnification of the image through the computer camera software, at this time the microscope only Power button and Light wheel available.

③**Charge Mode:** When the microscope is low on power, you can use the microscope while charging, use Type-C data cable to connect the microscope and the computer side, the microscope is turned on, select Charge Mode, you can enter the charging state. (Note: you can also use daily cell phone charger, computer, rechargeable battery to charge the microscope)

Work with TV


Operating steps

Step 1. Connect this microscope to the TV via HDTV cable

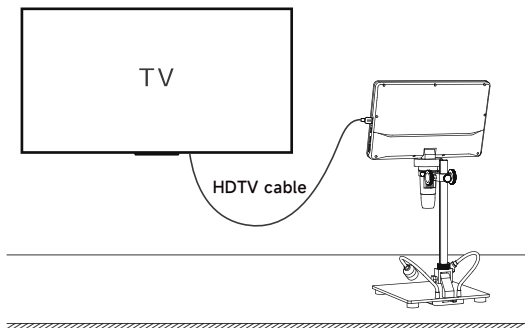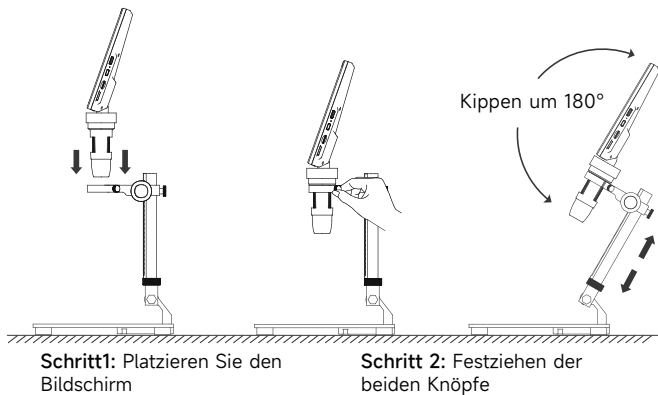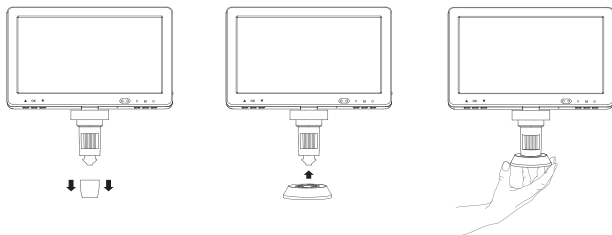
Step 2. Turn on the microscope and the image will appear on the display. At this time, microscope button and remote control button can be operated; when operating photo / video, the file will be automatically saved in TF card.

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